

August 2009

Advantage  
Business Media

www.surgicalproductsmag.com


# SURGICAL PRODUCTS

Your Leading Source for Cutting Edge Surgical Technology

## 2010 BUYERS' GUIDE

*The Complete Product  
& Supplier Directory for  
the Surgical Market*

- *Equipment, Systems,  
Services and Supplies*
- *Directory of Suppliers*

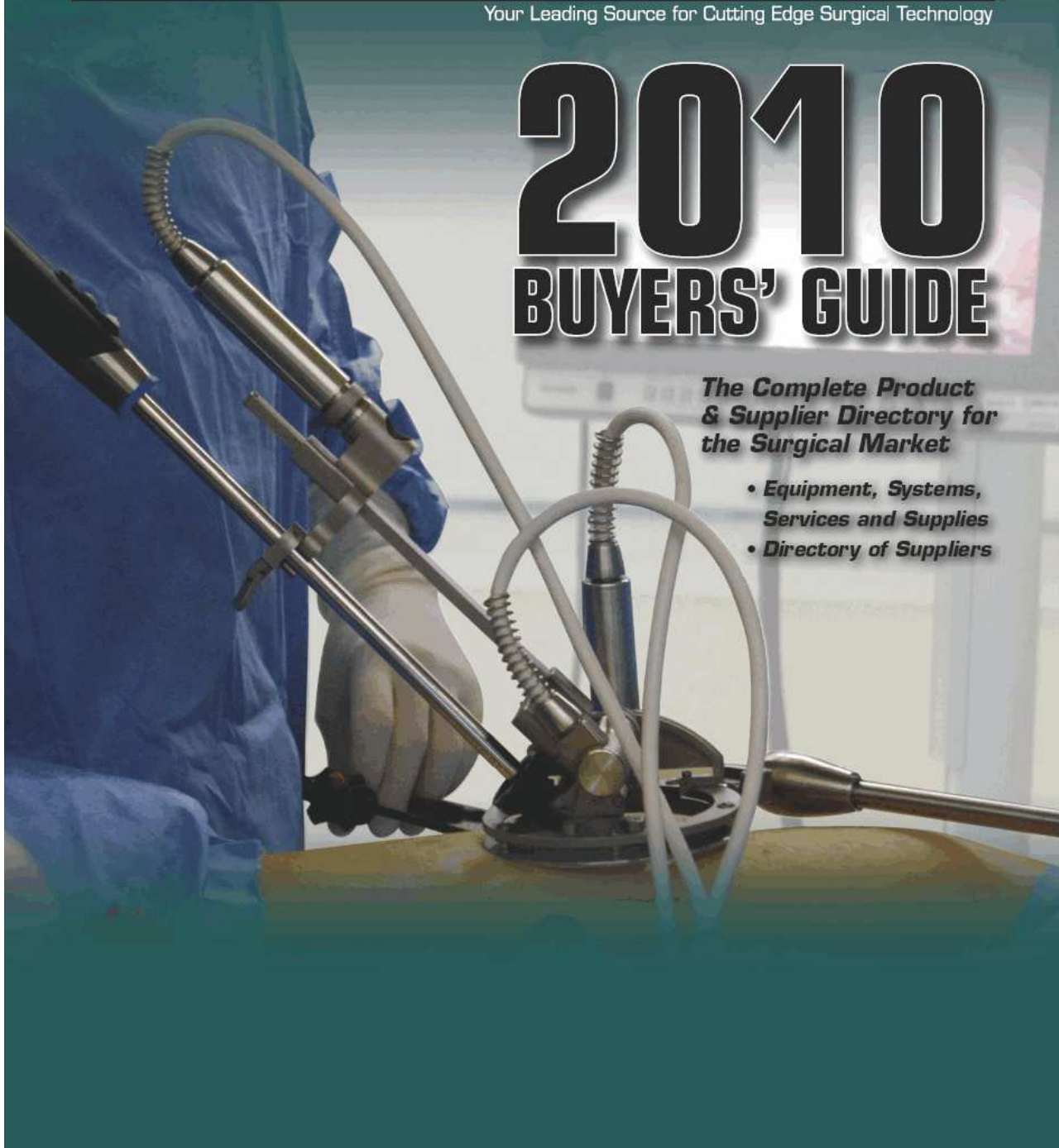



## The most comprehensive source for STERIS Surgical Table Accessories



# SURGICAL PRODUCTS®

### Now On Surgical Products

- General Surgeons On List Of...
- Hospital Adopts UVC Lights...
- Procedure Trainer With Removable...

**FREE** SUBSCRIPTIONS  
Newsletters • Magazine

### Subscriptions

- Home
- News
- Features
- ESP Awards
- First Cuts
- Products**
- Editor's Choice
- Capital Equipment
- Instrumentation
- Infection Control
- Patient & Staff Safety
- Wound Closure & Healing
- Magazine**
- Digital Edition
- Subscribe
- Resources**
- Calendar of Events
- Buyer's Guide
- White Papers
- Free Buyer's Guide Form
- Media Kit
- Contact**
- Staff
- Feedback

### Editor's Choice

#### Hands-Free Vision

The ViKY® System "Vision Kontrol for endoscopy" from **EndoControl** is a motorized endoscope holder for laparoscopic and thoracic surgeries. The holder features:

- Ease of installation with positioning just above the patient's abdomen and held by a passive arm.
- The ability to be fixed anywhere along the rail of the operating table.
- The ability to be moved while the system is in place.
- 3 motors of the camera driver allowing 3 degrees of freedom for the endoscope (2 rotations to visualize the entire abdomen cavity and an axial translation for zooming).
- Control by a foot switch command.
- The ability to be freely moved by hand in any direction (right/left, up/down and zoom) when disabled.
- The ability for the endoscope to be easily removed anytime for cleaning.
- Full immersability of the camera driver and passive arm with the ability to be sterilized by autoclave.



For more information, visit [www.endocontrol-medical.com](http://www.endocontrol-medical.com)

*Surgical Products* is conducting a brief, one-question survey about healthcare insurance coverage. [Click here to answer the question.](#) Thank you!

#### EndoControl

5, av du Gd Sabion, 38700la Trousche

- [Email for more information](#)
- [Company's other products](#)
- [Email to a colleague](#)
- [See similar products](#)
- [Printer friendly format](#)



© 2009 Advantage Business Media. All rights reserved.  
Use of this website is subject to its terms of use  
[Privacy Statement](#)